



Eaze **WAVE**

OPERATING GUIDE

Manufactured in the USA

Table of Contents

How EAZE WAVE Works	02
What is Included with the System	04
EAZE WAVE Pad Set Up	05
EAZE WAVE Controller Set Up	08
EAZE WAVE Controller Output Settings Guide	10
Adjusting Output Settings	11
EAZE WAVE Operating Instructions	12
Specifications and Disclaimers	13
Troubleshooting	15
Care Instructions	16

How EAZE WAVE Works




The EAZE WAVE System implements a patented stimulation protocol to provide relaxation prior to and during sleep.

The frequencies used are all within the audio spectrum, but the non-audible EAZE WAVE technology creates a noise-free environment where these relaxing frequencies are converted into low-intensity magnetic waves.

Stimulating the Body's Energy Points


The EAZE WAVE system comprises of a EAZE WAVE Controller that delivers a narrow band radio frequency signal at 28 Kilo Hertz with side bands that incorporate relaxation frequencies into a conductor that is embedded within an EAZE WAVE Pad. This radio frequency signal generates low intensity magnetic waves as it transverses down the conductor that stimulate key energy points along the spine positioned at L4, Th10, Th3, C6 vertebrae locations. The aim of this stimulation is to balance the body's energy and improve the flow of energy along the spine to promote relaxation prior to and during sleep.

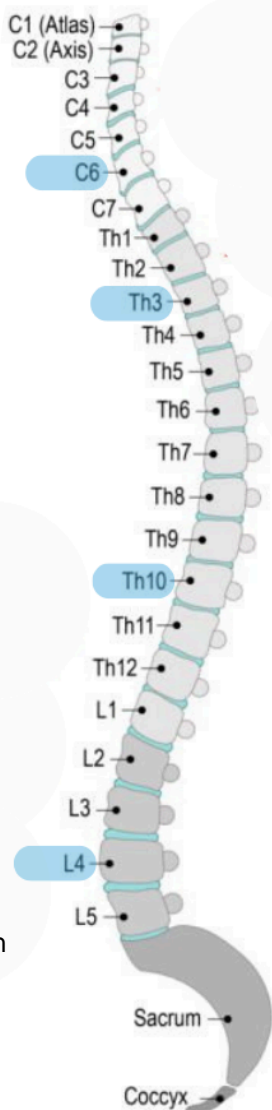


Spinal Model with Vertebrae Location

What Is Included with the System


Inside the EAZE WAVE box, you'll find everything you need to keep your system running smoothly:

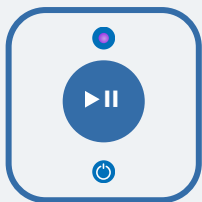


**EAZE WAVE
Controller**

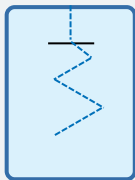

**EAZE WAVE
Fitted or
Travel Pad**

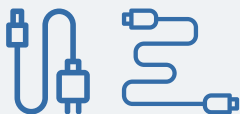

**Power Cord
and
Controller
Cable**

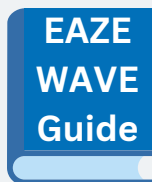

**EAZE WAVE
Operating
Guide**

EAZE WAVE

Pad Set Up

Follow the steps below to correctly set up EAZE WAVE Fitted or Travel Pad on your mattress

1


Unfold
EAZE WAVE Fitted or Travel Pad

2

Connect the controller cable to the
Fitted or Travel Pad connector

Image 1

Image 1

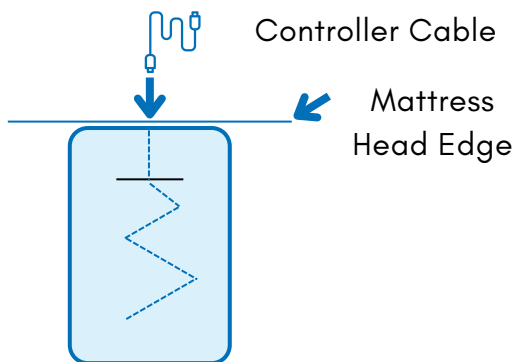

3

Position the Fitted or Travel Pad flat directly onto the mattress and underneath the mattress protector. Ensure the pad's connector end is on the mattress head end, not the foot end.

4

The Fitted Pad will self-align. For the Travel Pad, align the Pad along the side and top edge of the mattress so that the connector end overhangs the mattress head edge.

5

When positioning the Fitted or Travel Pad on the mattress, ensure that the black line is visible, which is used for positioning the pillow on the bed.

6

Place all bed items (mattress protector, fitted sheet, etc.) over the Fitted or Travel Pad ensuring the Pad is flat on your mattress.



EAZE WAVE


Fitted Pad

Set Up

7

Place your pillow on the mattress, ensuring it sits in between the black line and the top edge of your mattress
(Image 2)

Image 2

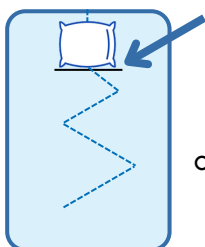

Black Line Pillow Alignment

Align the bottom center of your pillow with the black line

8

You are ready to move onto the next step:
EAZE WAVE Controller Setup
(next page)

EAZE WAVE Controller Setup

Follow the steps below to correctly set up and operate the EAZE WAVE Controller


1

Check to make sure the Power Button is Off (Button Out)

2

Connect the controller cable to the Controller (**Image 3**)

Image 3

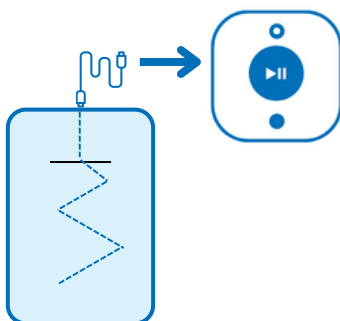


3

Connect one end of the Power Cord to the Controller, and the other end to a power outlet (**Image 4**)

4

Place the Controller near to your bed, on a flat surface, such as a nightstand

Image 4

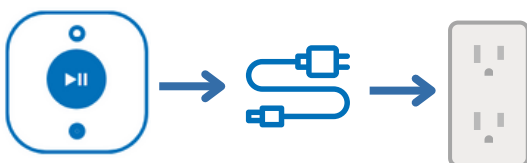

EAZE WAVE Controller Setup

5

Press Power On Button to turn the System On. You should see a **Blue Light** appear

6

Press the Activate Button and the light will change to **Purple**, indicating the system is now activated

EAZE WAVE Light

Blue Light: System is ready to use

Purple Light: System in activated

Activate/Pause

Output
Setting
Switch

Power On (Button In)
Power Off (Button Out)

EAZE WAVE

Controller Output Settings Guide

EAZE WAVE has four (4) output settings. See details below on how to use and adjust the output settings:

Setting 1

We recommend starting with **Setting 1** during the **first 60-nights** to allow the body to get used to the system.

Setting 2

After 60 consecutive nights, change to **Setting 2** if relaxation is not being achieved at **Setting 1**.

Setting 3

After an additional 90 consecutive nights, change to **Setting 3** if relaxation is not being achieved at **Setting 2**.

Setting 4

If relaxation is not being achieved at **Setting 3**, change to **Setting 4**.

Move at Your Own Pace

You may stay at a setting for as long as it helps you relax, even if it's longer than the recommended setting period. You can also move to a higher setting at a faster pace than recommended.


Adjusting Output Setting

The output setting switch is located on the bottom of the Controller.

SAFETY FIRST: To change the setting, first disconnect the power cord from the power supply and from the Controller to ensure that there is no power flowing through.

Then, flip the Controller upside down to access the switch panel.

Using the tip of a pen or similar object move the output setting switch (micro switch) to the output setting you desire.

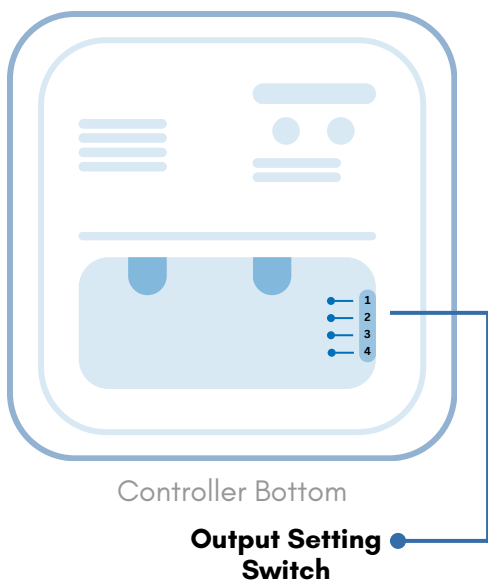

EAZE WAVE Operating Instructions

After setting up the EAZE WAVE Fitted or Travel Pad (pages 5 to 7) and the EAZE WAVE Controller (pages 8 to 11) you are ready to use the EAZE WAVE system.

Ready to Sleep

When ready to sleep (lights off, TV off), turn the Controller on by pressing the power button. When the controller is on the power button will be IN (page 9) and a **Blue** light will appear. Press the Activate button, and a **Purple** light will appear (page 9). Once the **Purple** light is on you may start experiencing the EAZE WAVE relaxation technology.

On Awakening

Turn the Controller off by pressing the power button. When the controller is off the power button will be OUT (page 9) and the light will also turn off.

Important

If you forget to turn off the Controller on awakening and the **Purple** light is on then before using EAZE WAVE again always power down the Controller and ensure that the light turns off before powering back up again.

When to Pause EAZE WAVE

If you need to leave your bed for more than 10 minutes while using EAZE WAVE, please pause it. Just press the pause button (page 9) when you get up and press the activate button when you return to bed. Ensure the system light is **Purple** when lying back down, indicating EAZE WAVE is ready to use.

Specifications and Disclaimers

Controller Specifications

V_{in}: 12V_{DC} regulated power supply

Bandwidth: 37Hz to 28KHz

V_{out}: 4.039V_{AC}

Carrier Freq: 28.4 kHz

I_{out}: 25.04μA

Disclosure and Disclaimer:

- Proper and standard use as provided in the OPERATING GUIDE is important.
- EAZE WAVE has not been evaluated by the Food and Drug Administration. This product is not intended to diagnose, treat, cure or prevent any disease.

Disclosure and Disclaimer:

- EAZE WAVE has been rigorously tested for safety. Specifically, compliant to IEC 60601-1 Medical Electrical Equipment (Part 1) General requirements for basic safety and essential performance and compliant to IEC 62233:2005- Measurement of electromagnetic fields of household appliances and similar apparatus with regard to human exposure.
- EAZE WAVE complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:
 - (1) This device may not cause harmful interference.
 - (2) This device must accept any interference received, including interference that may cause undesired operation.

Recommendation

To help your body adjust to the system, make sure to stay hydrated and limit the intake of alcohol and caffeine before going to sleep.

TROUBLESHOOTING

Identifying Possible Issues and Troubleshooting:

- On awakening, the Controller light should be **Purple** (unless the system was paused). A **Blue** light on awakening indicates a power drop or surge during the night. Restart the system.
- A flashing **Purple** light indicates an issue reading the frequencies.
- The issues above are **rare** and should fix themselves by restarting the system.
- If restarting the system doesn't fix these issues, please contact us at **info@eazewave.com**

How to restart/reset the System?

Turn the system Off by pressing the Power Off Button and turn the system back On by pressing the Power On Button.

CARE INSTRUCTIONS

**Due to the Electronic Nature of the
EAZE WAVE Pad**

DO NOT WASH
DO NOT DRY-CLEAN

SPOT CLEAN ONLY

by blotting the stain. Apply a stain remover to a clean white cloth and blot the stained area. Use another clean white cloth dampened with water and blot again to rinse. Let the cover dry completely. Repeat the process if needed.



EAZE WAVE
San Diego, California
eazewave.com
info@eazewave.com